Gebrauchsanweisung  
Instructions for Use  
Mode d'emploi

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umpe

akkubetrieb

► = Info

➔ = Eingabe

PEGA® COMBI LINE Y (REF 10523Y)

**PEGA® COMBI LINE Y**

CE 0123

**VENNER**

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**This page left intentionally blank**

**Intended Use**

PEGA® COMBI LINE Y consists of an infusion tube with a sterile filter and a peristaltic tube for parenteral administration of drugs in ambulatory and in-patient infusion therapies.

PEGA® COMBI LINE Y consists of a special check valve that prevents back-and freeflow. A second connector with siphon (Y-connector) allows additional medications to be introduced.

PEGA® COMBI LINE Y can be used with PEGA® BAG, or it can be used with PEGA® SPIKE as a connection between PEGA® COMBI LINE Y and infusion bags with pierce-able septum.

PEGA® COMBI LINE Y is only intended for use with PEGA® Infuse, PEGA® PCA, PEGA® plus and PEGA® light infusion pumps.

**Indication for Use**

The PEGA® COMBI LINE Y is suitable for intravenous, subcutaneous, epidural, and

intrathecal administration of drugs from a PEGA® infusion pump.

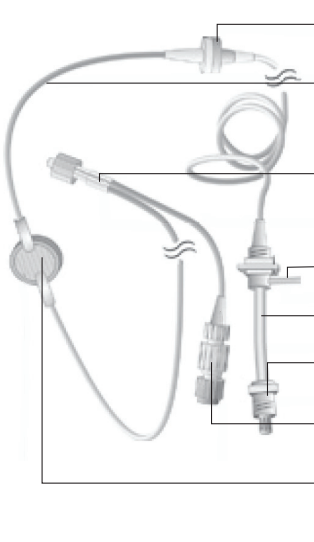
View	No.	Description / Specification
	1	Siphon / Check valve
	2	Infusion tube      L = 2,070 mm, ID = 0.5 mm, Dead volume PEGA® COMBI LINE Y V <sub>0</sub> = 1,600 µl
	3	Connection (patient) / Luer-Lock (male)
	4	Segment of the pump (Silicone) consisting of Pressure tube (4) and Peristaltic tube (5)
	5	
	6	Connection (reservoir) / Luer-Lock (male)
	7	Y-connector Luer-Lock (female) with siphon
	8	0.2 µm Sterile filter 0.1µm PVDF air-out-filter with embossed device LOT number

Fig. A: View PEGA® COMBI LINE Y

## English

### Contraindication

Do not use drugs which are not compatible with the material (PVC, PES, PVDF, PMMA, ABS, silicone rubber) of the PEGA® COMBI LINE Y.

### General Safety Instructions

- ◆ Please consult the manual of your PEGA® pump before using a PEGA® COMBI LINE Y.
- ◆ This disposable item is intended to be used on a patient only once in order to guarantee optimum performance and hygiene.  
  
If the item is used again, even after reconditioning, any liability on the part of Venner Medical (Deutschland) GmbH shall cease to exist.
- ◆ Please also respect the manuals of the PEGA® BAG and the PEGA® SPIKE.
- ◆ Carry out all of these steps on a clean and low germ concentration surface.
- ◆ Hands must be disinfected previously. Do not touch open luer-locks (danger of contamination).
- ◆ Use only original PEGA® accessories.
- ◆ In the case of parallel infusions, be sure to connect an additional check valve at the patient end of the infusion tube to prevent different medications mixing in the tube.

### Additional Information

PEGA® COMBI LINE Y is latex, DEHP and pyrogen free.

- ◆ PEGA® COMBI LINEs Y are provided in sterile packages. Do not use any disposable if the package appears damaged or opened.
- ◆ Use PEGA® COMBI LINEs Y only once and discard the disposable as per local regulations.
- ◆ Only remove the vented cap of the PEGA® COMBI LINE Y when connecting to the PEGA® BAG / PEGA® SPIKE.
- ◆ PEGA® COMBI LINE Y has to be substituted with every change of reservoir and must not be used longer than 30 days.
- ◆ A maximum delivery volume of 500 ml must not be exceeded.
- ◆ The pump with the PEGA® COMBI LINE Y must not be positioned higher than 1 m above the heart of the patient during the therapy.
- ◆ The PEGA® pumps are correctly closed only if the clips are fully engaged. If the clips are not correctly engaged a safe operation of the infusion pumps cannot be warranted.

### Confirmation of sterilization

Yellow Color indicator for EO sterilization method on the transparent side of the sterile

pouch advises that the sterilization is done and passed.

## Inserting

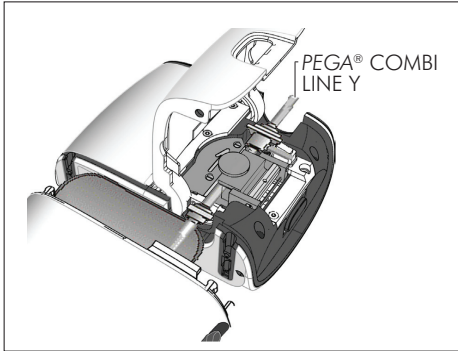
Fig. B: view with infusion pump (example PEGA® Infuse)

1. Open the rear pump cover.
2. Open the locking clip.
3. Open the sterile packaging of the PEGA® COMBI LINE Y.
4. Insert the peristaltic tube of the PEGA® COMBI LINE Y into the upper bearing of the pump. Ensure that the

## Priming

Before connecting it to the patient, prime the catheter system following the instructions in the manual of your PEGA® pump.

Ensure that no air bubbles are left in the tubing or the filter on the patient side (ribbed side). If this is not the case, flick the filter and tubing lightly with your finger until all the air bubbles are released.

### **Warning!**

The product is designed for a maximum pressure of 1,800 hPa!

5. Pull carefully on the Luer-lock-connection of the peristaltic tube so that this engages into the lower bearing of the pump.
6. If necessary, press the tube downwards with your finger so that it rests deep in the shaft.
7. Close the locking clip. The pressing block presses the peristaltic tube into the peristaltic drive.
8. Place the pump cover onto the pump and close it. The clips of the housing (one or two clips, depends on the PEGA® Pump version) must be heard to clip into place.

Continue the priming procedure until a drop of medication has formed at the end of the catheter system.

### **Caution**

Please note that the second connector is not primed automatically. It is essential that you prime the connector before use and also use a filter (not supplied).

## English

### Purchase Order

Purchase only complete boxes

REF.	Signification	Units per Box
10523Y	PEGA® COMBI LINE Y	10

Table 1: purchase order

### Conformity

Quality and safety of our products are the highest priority at Venner Medical. If a complaint should occur, please contact your responsible medical device consultant or

directly Venner Medical. In case of incidents with the PEGA products, report them to the competent authority in your country (e.g. BfArM in Germany).

**Meaning of the symbols on product and transport packaging**

 <b>MD</b>	Medical Device	 <b>REF</b>	Article number
	Use by date		Open here
	Single use only. Do not re-use		Keep dry
	Do not re-sterilize		Temperature limit (+5° to +40°C)
	Do not use if package is damaged		Keep away from sunlight
	Caution or Warning		Does not contain phthalates (DEHP)
	Caution, consult instructions for use		Does not contain or no presence of natural rubber latex
	Date of manufacture and country of manufacture		Pyrogen-free
	Batch code		Manufacturer
	Sterilisation with ethylene oxide		
	Single sterile barrier system		
	CE mark with code of notified body - according to MDD 93/42 EEC and MDR		

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